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
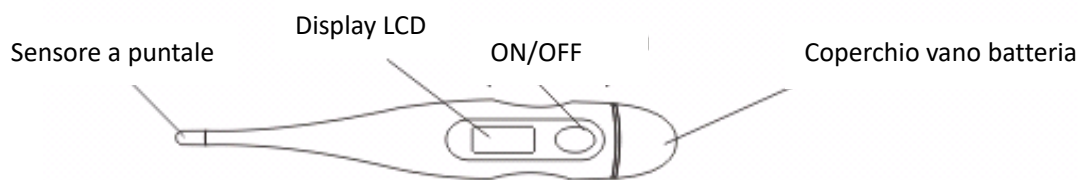
digiT-40 (rigid tip) ed digiT-40F (flexible tip) DIGITAL THERMOMETER

INTENDED USE:

Device **digiT-40 / digiT-40F** measures body temperature, which is measured at the metal measuring tip. This digital thermometer can be used under the tongue (oral), in the armpit (axillary) and in the rectum (rectal) to take temperature measuring results. An acoustic signal in "Bi, Bi, Bi..." will beep out once the temperature measuring process is complete. To ensure the best results and long-term satisfaction with this Digital Thermometer, we recommend that you read the instruction carefully before use and every time you will need in the future

SCHEMATIC DIAGRAM:

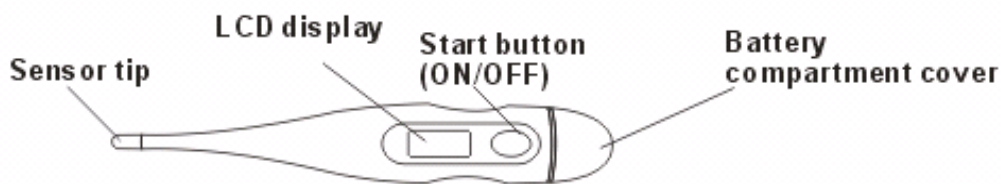NOTICE:

- During the measurement, do not walk, run or talk.
- The thermometer should only be used under the supervision of an adult.
- Do not store the thermometer where it will be exposed to direct sunlight, dust or humidity. Avoid extreme temperature.
- Dropping or subjecting your thermometer to strong shocks should be avoided.
- Do not disassemble the thermometer, except to replace the battery.

CAUTION:

- The thermometer contains small parts (e.g. battery), which could be swallowed by children. Therefore, keep the thermometer inaccessible for children.
- If the device will not be used for an extended period of time, make sure you remove the battery when you store it.


OPERATING INSTRUCTION

RECTAL (WAY) MEASUREMENT:

This is the most accurate method from a medical point of view, because it comes closest to the core body temperature. The thermometer tip is inserted carefully into the rectum for a maximum of 2 cm, after having previously applied a water-soluble lubricating solution to the orifice and the tip of the thermometer. Never force the insertion of the thermometer in case of resistance. Keep the thermometer in position until the end measurement beep alarm is emitted. **The usual measuring time is approximately 40 seconds.**

ORAL (WAY) MEASUREMENT:




There are different heat zones in the mouth. To ensure that reading is as accurate as possible, place the thermometer tip to the left or right of the root of the tongue. The thermometer tip must have constant contact with the tissue during the reading and be placed under the tongue in one of the two heat pockets at the back. Keep the mouth closed during the reading and breathe evenly through the nose.

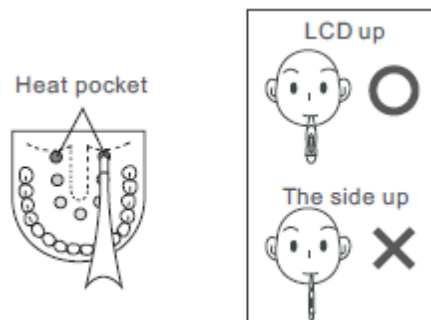Do not eat or drink anything before the measurement (30 minutes before do the measurement at least). **The usual measuring time is approximately 60 seconds.**

AXILLARY MEASUREMENT:

Placing the thermometer in the armpit provides a measurement of surface temperature. The usual measuring time for this method is **approximately 90 seconds**. It should be noted, however, that an exact reading cannot be obtained if, for example, the armpits have been allowed to cool.

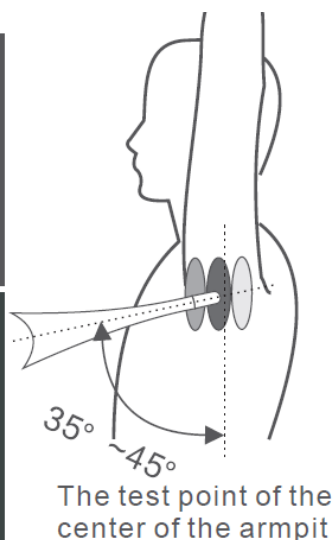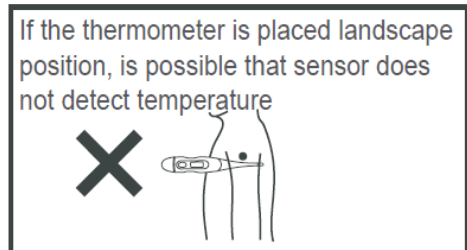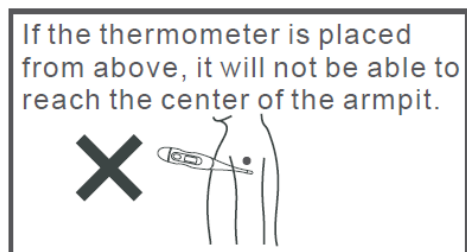Under this case, we recommend extending the measuring time by around 5 minutes in order to obtain the most precise possible reading that corresponds as closely as possible to the core body temperature (this independently by end measurement beep sounds).

Note: We strongly recommend the rectal method as the most accurate method for identifying the basal temperature, and advise you to extend the measuring time by 3 minutes after the beep (independently by measurement mode).

OPERATING:

- Press the ON / OFF button to turn on the thermometer. After a short beep indicates that the thermometer is now switched on. Simultaneously, a test display is performed in which all the symbols

to illuminate the display.

- The thermometer is ready to use when the "Lo" and a flashing "° C" appears in the display.

Notes: In an ambient temperature below 32 ° C will show "Lo ° C," appears in an ambient temperature above 42.9 ° C "Hi ° C".

During the measurement, the current temperature will be displayed continuously. During this period, the symbol "° C" blinks. The measurement is ended when a constant temperature is reached. The end of the measurement signal 10 beeps. The highest measured temperature is displayed.

To extend battery life, turn off the thermometer after use by pressing the ON / OFF button again. Otherwise the thermometer will automatically turn off after approximately 8 minutes.

MEMORY FUNCTION:

- Press the ON / OFF button to turn on the thermometer. The last measured temperature is automatically saved is displayed. The temperature appears with a small "M" in the display.
- You can now perform another measurement or turn off the power with the on / off button.
- The stored value is automatically overwritten when a new measurement is started.

CLEANING AND DISINFECTION:

To prevent infections the thermometer should be cleaned thoroughly after each use. We therefore recommend the use of commercially available disinfectants (e.g. a mild solution containing alcohol).

BODY TEMPERATURE:

Normal body temperature varies within a certain range. The following table shows the normal range depends on the measurement method. Therefore, measurement results can be directly compared different measuring points either.

Normal temperature ranges in ° C:	
Axilla	34.7 - 37.3

Oral	35.5 - 37.5
Rectal	36.6 - 37.7

However, the normal range of a persons subject toother slight variations. Therefore, it is important to determine the respective normal range for your self. Run the temperature measurement by your self andother family members from - ideally in a healthy state - to determine the normal range for each individual.

BATTERY REPLACEMENT:

Replace the battery when the battery level indicator  shows low battery status.

USE ONLY BUTTON BATTERIES 1.5 V DC, type LR41 or AG3


1. Place the thermometer with face back up and pull-out the battery cover by pressing it slightly downwards and at the same time pulling it outwards. CAUTION, during this operation the battery could fall to the ground easily if you apply excessive energy.

2. Take out the battery holder for approx. 18mm, take out the battery by a tool and replace it with a new one

3. Put back the battery cover in the right direction and make sure that it is properly locked in place.

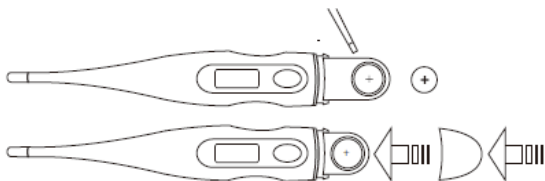INSTRUCTIONS FOR MAINTENANCE:

This device don't need of particular procedures of maintenance during its expected life except that battery replacement. Life expectation for this device is 3 years from the manufacturing data or 2 years from first use data (without batteries contained and in respect of storage conditions limits). After this period, is recommended to use device after an accuracy check done from qualified third part institute.

**SPECIFICATIONS:**

Model No.:	digit-40 (rigid tip model)// digit-40F (flexible tip model)
Display:	LCD,3.5digits, 16mmx7.5mm
Temperature measurement:	thermistor
Apply for:	Rectal / Oral / Axillary
Temperature measuring time:	40/60/90 seconds
Measuring range:	32.0 °C ~ 42.9 °C
Accuracy:	+/- 0.1 °C from 35.5 °C to 42 °C +/- 0.2 °C in the rest measuring range
Resolution:	0.1 °C
Memory records:	1 set
Low voltage Indication	YES
Automatically turn-off function	YES (after 8 minutes)
Power source and voltage:	DC 1.5V button batteries, type LR41 or AG3
Battery life:	Approx. 1200 times or 2 years
Dimension:	126mm×20mm×11mm
Unit weight:	Approx. 9.4g (digit-40) – 9.2g (digit-40F) including battery
Work environment:	Temperature 5°C to 40°C; Humidity: 15%RH~85%RH; Atmospheric pressure: 70KPa~106KPa
Transport/Storage condition:	Temperature -20°C~55°C; Humidity: 10%RH~85%RH; Atmospheric pressure: 70KPa~106KPa

*If there is any specification change for the product, these will applied without prior notice.

SYMBOL EXPLANATION:

1		Caution
2		Keep dry
3		Read the instruction before use
4		Type BF Applied Part Device
5	IP54	Degree of protection of an electrical device against accidental or intentional penetration of foreign solid bodies, and protection against the penetration of liquids. 1st Number - Protected against dust (no harmful deposit) 2° Number - Protected against water sprayed from all directions
6		Low battery status
7		Date of Manufacturer
8		Manufacturer

CA-MI S.r.l. - Via Ugo La Malfa, 13 – Frazione Pilastro, 43013 Langhirano (PR) – Italy

This product is MADE IN CHINA

THE THERMOTER WAS DESIGNED AND MANUFACTURED IN COMPLIANCE WITH THE FOLLOWING EUROPEAN STANDARDS:

IEC 60601-1	Medical electrical equipment – Part 1: General requirements for basic safety and essential performance
IEC 60601-1-2	Medical electrical equipment – Part 1: General requirements for basic safety and essential performance – Collateral standard: Electromagnetic compatibility – Requirements and tests
ISO 80601-2-56	medical electrical equipment part 2-56: particular requirements for basic safety and essential performance of clinical thermometers for body temperature measurement
ISO15223-1	Medical devices--symbols to be used with medical device labels, labelling and information to be supplied lied--part1:general requirements

This product complies with the European Directive for Medical Devices 93/42/EEC and carries the CE mark.

Electromagnetic compatibility:

The device complies with the IEC 60601-1-2 standard for electromagnetic compatibility. Please do not let thermometer become source of the disturbance, also try your best avoid source of the disturbance. Inquire at the Manufacturer for details on this measurement data. We reserve the right to make technical and design changes in the course of continuous product improvement .

WARNINGS!

- This device should not be used in the vicinity or on the top of other electronic equipment such as cell phone, transceiver or radio control products. If you have to do so, the device should be observed to verify normal operation.
- The use of accessories and power cord other than those specified, with the exception of cables sold by the manufacturer of the equipment or system as replacement parts for internal components, may result in increased emissions or decreased immunity of the equipment or system.

Warranty

This product is covered by warranty for 2 years from the date of purchasing. Damage resulting from incorrect use or abuse is not covered by the warranty. Battery and packaging are excluded from the warranty. If you find that the thermometer is defective and not operating correctly, please firstly check the battery before sending in for repair. Warranty is voided in case of tampering, if the defect is caused by improper use or if the damage is not due to manufacturing faults (accidental drop, inaccurate transportation, etc.). Warranty is valid from date of purchasing (certified from Invoice or receipt).

IMPORTANT INFORMATION FOR CORRECT DISPOSAL OF THE PRODUCT IN ACCORDANCE WITH EC DIRECTIVE 2012/19/UE-WEEE: The symbol on the device indicates the separated collection of electric and electronic equipment. At the end of life of the device, don't dispose it as mixed solid municipal waste, but dispose it referring to a specific collection center located in your area or returning it to the distributor, when buying a new device of the sample type to be used with the same functions. This procedure of separated collection of electric and electronic devices is carried out forecasting a European environmental policy aiming at safeguarding, protecting and improving environment quality, as well as avoiding potential effects on human health due to the presence of hazardous substances in such equipment or to an improper use of the same or of parts of the same. Caution: The wrong disposal of electric and electronic equipment may involve sanctions.

DISPOSAL OF SPENT BATTERIES (Directive 2006/66/EC): This symbol on the product indicates that the batteries should not be treated as household waste. Ensuring that the batteries are correctly disposed of contributes to preventing potential negative consequences for health and the environment, which could otherwise be caused by their inappropriate disposal. Recycling materials helps to preserve our natural resources. Take the spent batteries to collection points for recycling. For more information about collection and recycling of spent batteries or the product, please contact your Municipality, your local waste disposal service or the shop where you purchased the device.

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